

DIAPHRAGM PULSATION DAMPER

FPD 06, 1.06, 10, 1.10



FPD 06, 10 KP.51 Z



FPD 06, 1.06 Z



FPD 10, 1.10 Z

Concept

The KNF diaphragm pulsation damper can effectively reduce the pulsation after the pump by up to 97%, whereby it is important to note that the type of pulsation damper selected and the correct installation are essential in ensuring an optimal performance.

The pulsation damper is made up of only a few parts and is constructed in a simple but effective way. A diaphragm separates the liquid section (upper part) from the damping chambers (lower part). The damping elements (situated in the damping chambers) combined with the separating diaphragm remove a large part of the pulsation caused by the pump. In the upper part of the pulsation damper there are also two choke orifices which aid the damping process.

Features

High chemical resistance

The use of chemically resistant materials such as PP, PVDF, PTFE, FPM or other material combinations for the parts which come in contact with the liquid allows almost all neutral or corrosive liquids to be pumped.

Durable and maintenance free



The carefully considered design of this very compact pulsation damper ensures safe operation and a long life even under the most severe conditions.

Areas of use

The KNF diaphragm pulsation damper is very versatile thus allowing it to be used not only in conjunction with a KNF pump, but in many general situations where less pulsation is required.

The pulsation damper can be used in the following areas:

- To reduce the vibration in hoses and pipes.
- To help remove pulsation which is preventing the system from functioning correctly.
- The protection of instrumentation connected after the pump.

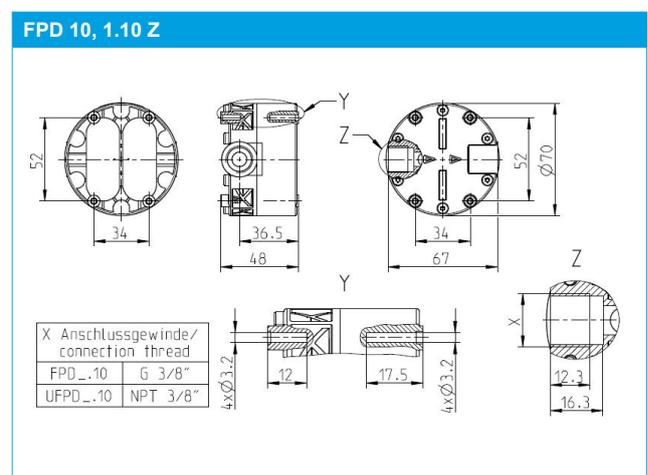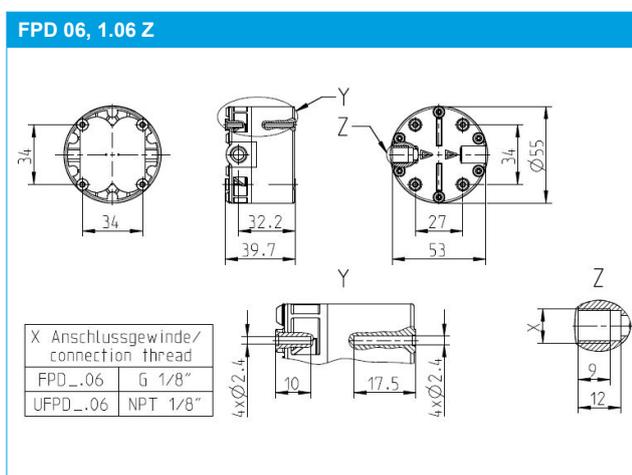


THE KNF MODULAR CONCEPT OF SELECTION

Our versatile self-selection program allows you to personally determine the optimum characteristics that you require from your pulsation damper. Select your diaphragm pulsation damper from the following characteristics:

KNF pump type being used	Max. allowed pressure (mWg)	Pulsation damper type			
		Basic model	1	2	accessory
NF 5, NFB 5, NFB 25, NF 10, NF 25, NF 30, NF 60, NF 100	20	FPD 06			Z
NF 1.5, NF 1.10, NF 1.30, NF 1.60, NF 1.25	60*	FPD 1.06			Z
NFB 100, NF 300, NF 600	20	FPD 10			Z
NF 1.100, NF 1.300, NF 1.600	60*	FPD 1.10			Z
	Example	FPD 06	KP	.51	Z

* For the KT- and TT version: max. allowed system pressure = 40 mWg

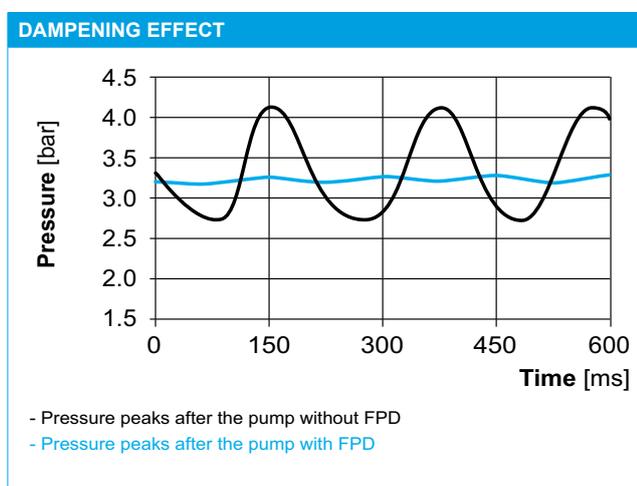

Selection and operation

- Hydraulic connections:
 - NPT thread
 - G thread
- The KNF diaphragm pulsation damper has been constructed so that the remaining pulsation will be lower than 200 mbar. This value can be improved upon if the system parameters (system pressure, motor speed, etc.) are known.
- There is the possibility that the pulsation damper will increase the pressure after the pump. This should be taken into consideration when choosing the pump.
- The pulsation damper should be located as close as possible to the pump outlet. Long connection hoses between the pump and the pulsation damper can create vibrations.

1 MATERIALS		
KP/ KP.51*	Housing Diaphragm	PP EPDM
KT	Housing Diaphragm	PP FFPM
TT	Housing Diaphragm	PVDF FFPM

* food grade conformity according to the standard NSF/ANSI 169, for further information see options.

2 OPTIONS	
.51	NSF certified / FDA conform



Options

NSF National Sanitary Foundation

NSF is market leader in the development and controlling of standards relevant for equipment used for handling food. By using different toxicological tests our products will be certified according to the standard NSF/ANSI 169. This certification will confirm that all of the products with the code .51 are certified for the use with foods/consumables. In addition to this it also means that the pumps are constructed with wetted parts which come with a FDA declaration of conformity. A yearly audit from NSF will be carried out to ensure that these standards are being maintained. A list of the various products are available on request.

* FDA = Food and Drug Administration